


HEINE QUALITY
MADE IN GERMANY

HEINE® iC2 Funduscope
Digital Ophthalmoscopy

NEW!

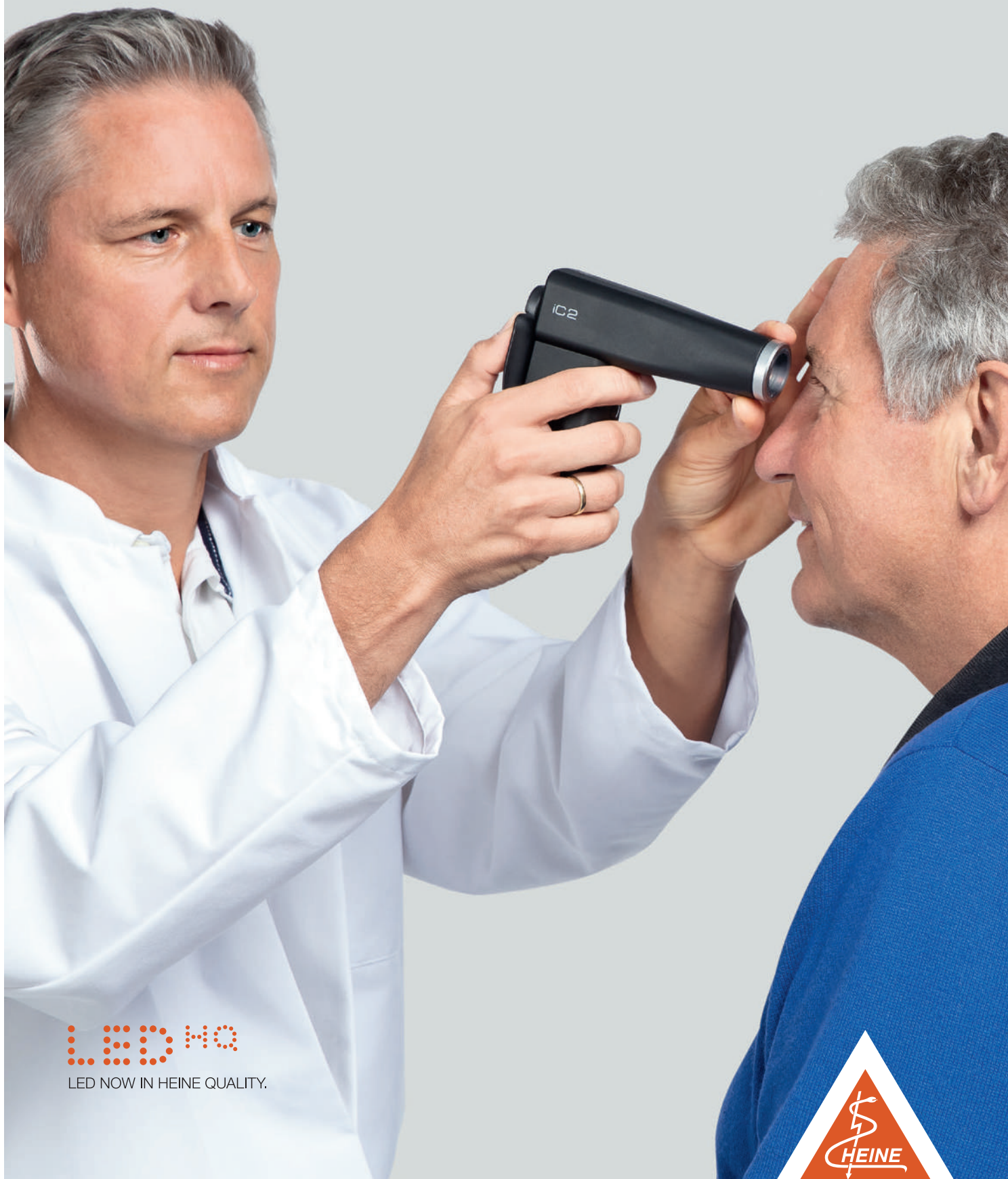

LED HQ

LED NOW IN HEINE QUALITY.



The new HEINE iC2 Funduscope is the next generation of direct, digital Ophthalmoscopes. Conveniently capture high quality fundus images.



✓ **Mobile and light**

Weighing only 300g, compact in size and supporting wireless operation, this instrument is ideal for convenient examinations even for patients who are confined to bed or have a disability.

✓ **Large field of view**

Field of view up to 34° (approx. 4x as large as the conventional direct ophthalmoscope with typically 9°).

✓ **Easy imaging of the periphery**

Scan the periphery as with a direct ophthalmoscope.

✓ **Sharp, shadow-free images**

Fine-tuned interaction with HEINE LED^{HQ} illumination with precision optics.



✓ **High-quality fundus images**

Early identification of diseases of the eye. Ideal also for consultations. Documentation of current condition and development over time. Excellent image quality with up to 12 MP**.

✓ **Compensation of patient's defective vision**


Stepless manual dioptric adjustment of +/-15 dpt. and auto focus with +/-3 dpt.

✓ **Photobiological Protection for the patient's eyes**


Complete brightness control and LED illumination which meets the highest international photobiological safety standards (DIN EN 15004-2, Group I)

✓ **iC2 App**

Simple mapping of images. E-mail and zoom functions. Storage of additional patient data. Always up to date thanks to updates from the App Store*.

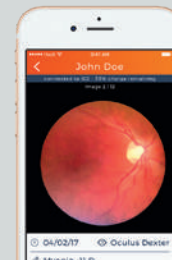


High-resolution images

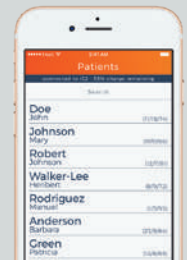

Select existing patients

Examples of fundus images with the HEINE iC2 Funduscope

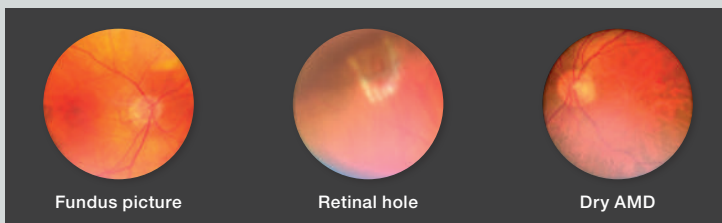

Explain findings to patients

Ordering information:

iC2 Set/8 for iPhone 8	
iC2 with mounting case smartphone iC 2/8 for Apple® iPhone® 8*, USB cord with E4-USB medical approved plug-in power supply and case	C-153.40.303
iC2 Set/7 for iPhone 7	
See above, but with mounting case smartphone iC 2/7 for Apple iPhone 7*	C-152.40.303
iC2 Set/6 for iPhone 6/6s	
See above, but with mounting case smartphone iC 2/6 for Apple iPhone 6/6s*	C-151.40.303
iC2 Set/5 for iPhone 5/5s and SE	
See above, but with mounting case smartphone iC 2/5 for Apple iPhone 5/5s and iPhone SE*	C-150.40.303

▶ Due to local regulations, some products may not be available in your country. For additional product configurations and further HEINE Quality Products please visit www.heine.com

* The Apple iPhone 8, 7, 6/6s, 5/5s or SE is not included in the scope of delivery. Apple, the Apple logo, iPhone and iPod are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

** With iPhone 6s, SE, 7 and 8



Имя HEINE характеризует продукцию высочайшего качества!

HEINE QUALITY
MADE IN GERMANY

"

-

"

<http://heine-med.ru/>
info@heine-med.ru

: (495) 902-59-26 .. , (495) 518-55-99

HEINE,

HEINE () -

!

